

PIEZON[®]

V-

CLASS

EMS⁺

ORIGINAL PIEZON® V-CLASS
NEW STANDARD



**NEW →
EMS SWISS
INSTRUMENT
V-PS**

> Suits nearly 90%
of all scaling applications



> Absolute precision
in the top 2 mm of
the instrument tip



> Conical shape
and short thread



> Ultra-slim
instrument design

**COPIES
WON'T DO IT**

- > Low prices do not make up for significantly shorter lifespan
- > Much faster wear and poorer treatment results
- > No EMS warranty when copies or non-EMS instruments are used

ORIGINAL PIEZON® V-CLASS
NEW STANDARD

NEW → ORIGINAL PIEZON® V-CLASS

EMS SWISS V-INSTRUMENTS^{PM}
AND V-HANDPIECES
FOR A NEW STANDARD
IN SCALING

- OPTIMUM TRANSMISSION OF POWER
- EASY HANDLING
- LONGER LIFE

**NEW STANDARD
OF EXCELLENCE**

ORIGINAL PIEZON® V-CLASS

V-CLASS DESIGN

V-CLASS INSTRUMENTS^{PM}

→ PROTECTED DESIGN

EMS SWISS V-INSTRUMENTS^{PM}
FOR ULTIMATE PRECISION, EFFICACY,
SPEED AND LONG LIFESPAN

QUICK AND SECURE

> Conical design for enhanced ergonomics. Instrument thread and sonotrode snap into place automatically. Half a turn is all it takes to secure the instrument – IT'S A SNAP


LONG LIFESPAN

> Conical design for reduced wear of instrument thread and longer lifespan in the practice – guaranteed water-tight coupling between instrument and handpiece


> Instrument made from a special alloy for improved lifespan

PRECISE AND EFFICIENT

> Conical design of instrument thread and sonotrode for perfect fit between instrument and handpiece – optimum transmission of power from piezoceramic drive to instrument tip

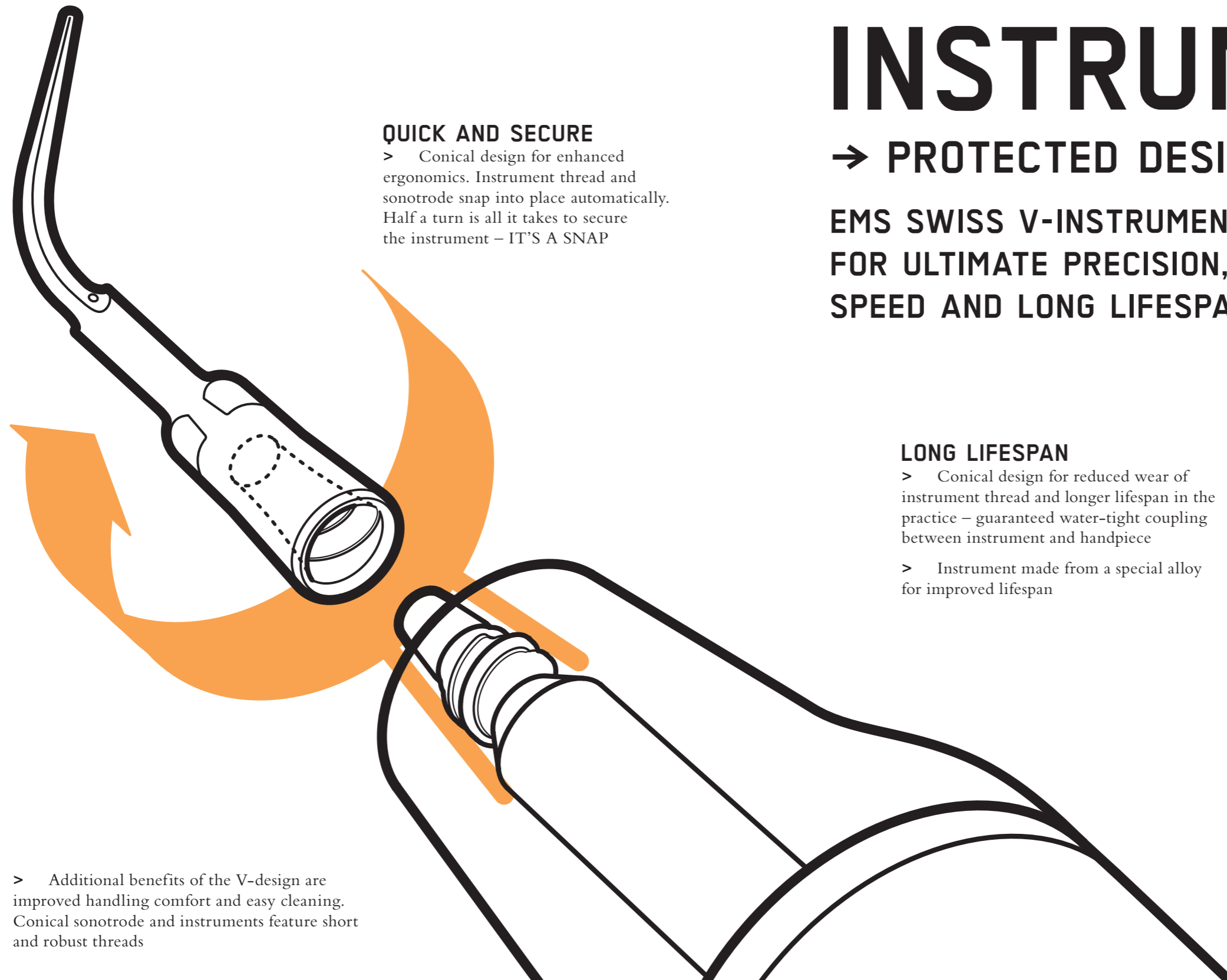

> Additional benefits of the V-design are improved handling comfort and easy cleaning. Conical sonotrode and instruments feature short and robust threads

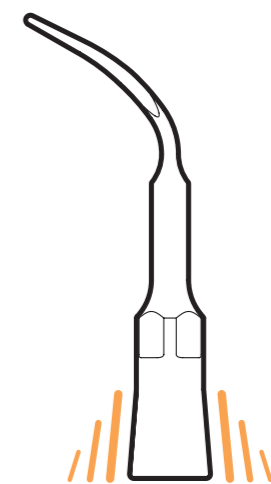

>



INSTRUMENT V-A

> For the removal of supragingival deposits in all quadrants



INSTRUMENT V-B

> For the removal of supragingival deposits in all quadrants



INSTRUMENT V-P

> For the removal of supra- and subgingival deposits in all quadrants



INSTRUMENT V-PS

> For the removal of deposits in all quadrants, even in deep periodontal pockets



INSTRUMENT V-PL1

> Left-curved instrument for periodontal debridement – especially suited for difficult access to approximal surfaces



INSTRUMENT V-PL2

> Right-curved instrument for periodontal debridement – especially suited for difficult access to approximal surfaces



INSTRUMENT V-PL3

> Instrument for cleaning and disinfecting periodontal pockets – especially suited for check-ups and recalls



> Download our brochure on unique NO PAIN technology at ems-swissquality.com



INSTRUMENT V-PL4

> Left-curved instrument with small ball [Ø 0.8 mm] on instrument tip for the periodontal debridement of furcations and concavities



INSTRUMENT V-PL5

> Right-curved instrument with small ball [Ø 0.8 mm] on instrument tip for the periodontal debridement of furcations and concavities



INSTRUMENT V-RT1

> Diamond-coated instrument for the entry location of root canals and removal of calcifications in the coronal third of the root system



INSTRUMENT V-RT2

> Diamond-coated instrument for the removal of shoulders and other obstacles, e.g. creation of straight canals for easy access to broken instruments



FILE HOLDER V-180°

> For use with RT3 instruments as well as Berutti and RE2



FILE HOLDER V-120°

> For the flexible use of different endodontic files, specially suited for anterior and premolar areas. Also suitable for RT3, RE2, Berutti and ESI instruments



FILE HOLDER V-90°

> For the flexible use of different endodontic files, including ESI instrument. Specially suited for the molar area

EMS SWISS V-INSTRUMENTS^{PM} → STANDARD WITH THE NEW-GENERATION AIR-FLOW MASTER PIEZON®, PIEZON® MASTER 700, miniMASTER® LED – OR UPGRADE

ORIGINAL PIEZON® METHOD → NEW GENERATION

NOW WITH EMS SWISS V-INSTRUMENTS^{PM}
AND V-HANDPIECES



AIR-FLOW MASTER PIEZON®

> Original AIR-FLOW® and Piezon®, all-in-one



PIEZON® MASTER 700

> Original Piezon® with NO PAIN technology



miniMaster® LED

> Original Piezon® Power Pack

PIEZON® V-CLASS UPGRADE → COMPLETE, SIMPLE

FOR ALL EXISTING AIR-FLOW MASTER PIEZON®,
PIEZON® MASTER 700 AND miniMASTER® LED

**V-CLASS
FOR ALL**



> EMS Swiss V-Instruments^{PM} A, P and PS including CombiTorque® in Steribox



> Original Piezon® V-handpiece LED



EMS-SWISSQUALITY.COM

EMS
ELECTRO MEDICAL SYSTEMS SA
Chemin de la Vuarpillièrre 31
CH-1260 Nyon

Tel. +41 22 99 44 700
Fax +41 22 99 44 701
welcome@ems-ch.com
www.ems-dent.com